

NEUATION



WITNESS
PERFECTION

iShak 3D-5 NXT

(Digital 3D Shaker)

iShak 3D-5 NXT (Digital 3D Shaker)

You strive for perfection and it's obvious that you need a device which offers the same level of reliable performance. Most existing 3D shakers in the market make tall promises of reliable performance and high capacity which they claim will make your work easier. But what happens, when you put that claim to the test? Made with low-quality internal components, these devices neither offer reliable performance nor a large loading capacity. Led by a weak motor, the performance is also lackluster leading to unsatisfactory user-experience.


Neuation's wide-ranging expertise in R&D helps shape our design and development decisions to better suit the user's requirements and expectations. We, also take the utmost care to use the highest quality raw materials and components in the manufacturing of our devices to ensure consistent, reliable and maintenance-free performance. iShak 3D-5 NXT, our Digital 3D Shaker comes with a maintenance-free brushless DC motor which delivers exceptional performance to complement the device's high loading capacity, resulting in thorough mixing of the solution even under big loads.



Product Features



Brushless DC Motor For
Quiet Performance And
A Maintenance-free Long Life.

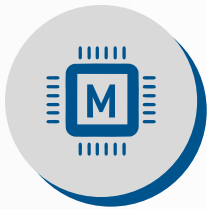

Microprocessor Controlled
Efficient Performance Even
Under High Loads.



5 Kg
Load Capacity



Variable Speed Setting
From 10 To 120 RPM

Variable Timer Setting From
1 To 999 Minutes Along
With An Infinite Mode.



Tumbling Angle Of 7°
Ensures Better Results



Pulse Mode Function
With Variable Timer And
Run-pause-run Function Ensures
Thorough Mixing Of Solution.



User-programmable
Mode Allows The User To
Set Upto 9 User-defined
Programs (Protocols).



Product Features



Digital Display



Non-slip Strength
Rubber Feet For
Maximum Stability.



Spiked Rubber Mat To
Ensure Stability Of The Tubes
While In Operation. Optional
Sticky Pad Mat Available.



Easy Access To The Device
Which Makes It Extremely
Easy To Clean And Maintain.

Robust design



IP21 Compliant Rating
Protects The Device Against
Penetration Of Liquids.



Safe Low
Voltage Device



Double Decker Design



Specifications

Specifications	iShak 3D-5 NXT
Weight Capacity	5 Kg
Variable Speed	10-120 RPM
Movement Type	Tumbling
Tumbling Angle	7°
Timer Range	1-999 mins & infinite mode
Operation Type	Timed & Continuous Operation
Pulse Mode	Yes (Run-Pause-Run Mode)
Plate Size	350 mm x 350 mm
Power Consumption	24W
Attachments	Spiked Rubber Mat (Included)
	Sticky Pad Mat (Optional)

Ordering Information

iShak 3D-5 NXT		
Product Name	Cat #	Details
iShak 3D-5 NXT	SH03-1200 EU/US/UK/AU	5 Kg 3D shaker Max RPM 120

Accessories	
Product Name	Cat #
Spiked Mat	SH03-M3-1
Sticky pad mat	SH03-M3-2

NEUATION

www.neuation.com

AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

*We also execute OEM projects. Specifications can be changed without notice for quality improvement.