

IT CAN'T GET BETTER THAN THIS



iSwix VT Digital Vortex Mixer

Vortex Mixers play a crucial part in any laboratory setup hence the need for reliable performance and durability, but most mixers have issues with lackluster performance owing to sub-par quality components and raw-materials used. Hence, these devices leave a lot to be desired and end-up being a waste of space and investment.

Developed after years of R&D, Neuation's best-in-class iSwix VT vortex mixer provides the best of performance due to its top-of-the-line internal components such as long lasting maintenance free BLDC motor. With class-leading features like inbuilt counter balance and 3-way switch, the device makes working even easier.



Product Features





Adjustable speed setting from 300 to 4200 RPM

Microprocessor based closed loop control for efficient performance even under high loads



Variable timer setting from 1 to 999 minutes and continuous mode





Orbital diameter of 4 mm

Pulse Mode function with Run-Pause-Run function for thorough mixing of solution





Digital Display shows the speed and time toggles after every 5 seconds

3-way switch for Touch, Continuous & Stand-by mode for unmatched convenience






Brushless DC motor for quiet performance and a maintenance-free long life

Product Features



Inbuilt counter balance ensures stable operation for complete peace-of-mind

Visual indication marker for max speed with attachments

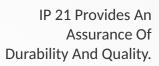









Small footprint saves valuable bench space and can also be used inside hoods and cold rooms

Standard Attachments



Rubber Disk

Attachment



Cup Attachment

Round Disc

Optional Attachments



Micro Centrifuge Tube

24 X 1.5/2ml







Flask

Specifications

iSwix VT	
Max Speed	4200 RPM
Orbital diameter	4 mm
Motor type	Brushless DC
Digital Display	Yes
Dimensions (W X D X H)	205 x 136 x 138.5 mm
Speed Setting	300-4200 RPM
Timer Setting	0-999 mins (min/sec display)
Pulse mode	Yes, Programmable Pulse
Max Load Capacity	500 gm
Weight	2.8 kg
Optional Attachments (not included)	Universal platform for attaching either - Microplate / Microtube / Flask / Round Disk Attachments

Ordering Information

iSwix VT		
Product Name	Cat #	Details
iSwix VT	VM-2110 EU/US/UK/AU	Digital Vortex Mixer
Attachment	VM-A1	Cup Attachment
	VM-A8	Rubber Cup Disk Attachment





www.neuation.com

AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

*We also execute OEM projects. Specifications can be changed without notice for quality improvement.