

ACCUMAX



A PERFECT FIT FOR SECURITY
HIGHLY-ADVANCED FOR PROTECTION

SAFETY CAP
Thread Size: GL 45

SAFETY CAP

Safety cap extracts solvents or mobile phases safely from the bottle or reservoirs without letting the hazardous vapour come out. It filters particulates in the air and prevents contaminants and dust entering solvents.



PRODUCT FEATURES



Advanced vapour valve keeps lab environment free from hazardous solvent vapour



Specially designed O-ring for air tight connection with glass bottle



To be used with solvents bottles for safer lab environment and lesser uncertainty of analysis



Chemically resistant against aggressive organic solvents

Efficient PTFE filter for preventing contaminants mediated noise




Freely rotating cap body prevents twisting tubes during solvent bottle exchange



High quality materials (PTFE and PP) ensure maximum safety in your laboratory



Blind fittings For closing unused ports on the cap

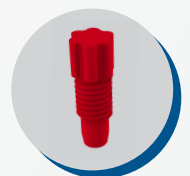

Ordering information:

HPLC Safety Cap			
Cat No	Product description	Pack size	Case size
SCPF-01	Safety Cap with filter valve 1 port & 1 filter	4 Cap, including 4 fitting & 4 plug	1
SCPF-02	Safety Cap with filter valve 2 ports & 1 filter	4 Cap, including 8 fitting & 8 plug	1
SCPF-03	Safety Cap with filter valve 3 ports & 1 filter	4 Cap, including 12 fitting & 12 plug	1
SCPF-04	Safety Cap with filter valve 4 ports & 1 filter	4 Cap, including 16 fitting & 16 plug	1
SCPF-05	Safety Cap with filter valve 5 ports & 1 filter	4 Cap, including 20 fitting & 20 plug	1
SCPF-A	Filter valve with 1 filter & 4 fitting & 4 plug	4	1

ACCUMAX

Accumax Lab Devices Pvt. Ltd.
www.accumaximum.com



ISO 9001:2015
ISO 13485:2016

*We also execute OEM projects. Specifications can be changed without notice for quality improvement.