


NEUATION

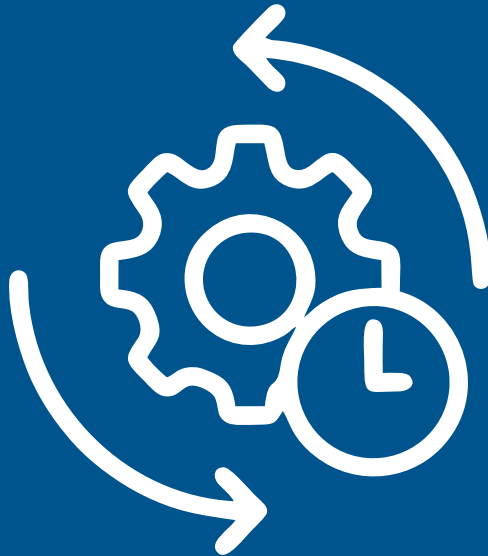

EFFICIENTLY FLEXIBLE
GEL CARD CENTRIFUGE

iFUGE GC12 Gel Card Centrifuge

iFUGE GC12 Gel Card Centrifuge

Gel card centrifuges are widely used in hospital and laboratory, for blood grouping test, ABO cross matching, blood routine inspection, Irregular Antibody screening, microcolumn gel, immunoassay and other tests.

The iFUGE GC12 Gel Card Centrifuge is a result of our extensive R&D into user's requirements and feedback. It is an efficiently flexible gel card centrifuge system that provides you freedom to select multiple suppliers of gel cards with replicating standard results. Considering the shortage of space in laboratories, the device is designed with a compact space saving table top footprint to occupy less space relative to its functionality. Class leading features and safety like Motorized latch, favourite mode, microprocessor controlled BLDC motor gives the user unprecedented convenience and ease-of-use.



Product Features




Easy to manage and program with minimal training requirement.



Capacity - 12 Gel cards



Digital Display with Speed & Time Setting

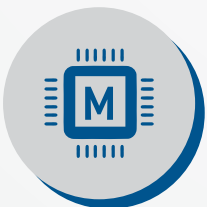

Microprocessor Controlled for Efficient Performance Even Under High Loads.



CE IVD Approved

Universal for all types of standard gel cards.



Max Speed- 2200 RPM



Timer Setting From 1 To 999 Mins & infinite



Lid Lock Safety Feature - Lid Opens Automatically On Run Completion



Product Features



One Touch RPM / RCF
Conversion & Short Spin
Operation Option on Keypad



Smart Airflow Design Ensures
Chamber Temperature Rise
Is Under Control



Favourite Mode:
Allows User to Set upto 99
User Defined Programs



Centrally placed latch
lock for effortless closure



Specifications

Specifications	iFuge GC12
Dimension (W x D x H)	381 x 480 x 185 mm
Weight of Equipment	8.3 kg (without rotor)
Permissible Temp.	5-40°C
Humidity	<80%
Altitude	Use upto an altitude of 2000 m above MSL
Pollution Degree	2
Environment	For indoor use only
Noise Level	≤ 65 dB
Power Adaptor Details	Input: 100-240 VAC, 50/60 Hz Output 24V $\overleftrightarrow{\bullet\bullet}$ 5A
Power Consumption	48 W

Ordering Information

iFuge GC12		
Product Name	Cat #	Details
iFUGE GC12	GC12 - EU/US/UK/AU	Gel Card Centrifuge 2200 RPM / 470g

NEUATION

*We also execute OEM projects. Specifications can be changed without notice for quality improvement.